



Rabbit Anti-D-DIMER antibody

SL3514R


Product Name:	D-DIMER
Chinese Name:	交联纤维蛋白二聚体抗体
Alias:	D DIMER.
Organism Species:	Rabbit
Clonality:	Polyclonal
React Species:	Human,Mouse,
Applications:	ELISA=1:500-1000IHC-P=1:400-800IHC-F=1:400-800IF=1:100-500 (Paraffin sections need antigen repair) not yet tested in other applications. optimal dilutions/concentrations should be determined by the end user.
Molecular weight:	340kDa
Form:	Lyophilized or Liquid
Concentration:	1mg/ml
immunogen:	homogenized human fibrin clot:Human homogenized fibrin clot.
Lsotype:	IgG
Purification:	affinity purified by Protein A
Storage Buffer:	0.01M TBS(pH7.4) with 1% BSA, 0.03% Proclin300 and 50% Glycerol.
Storage:	Store at -20 °C for one year. Avoid repeated freeze/thaw cycles. The lyophilized antibody is stable at room temperature for at least one month and for greater than a year when kept at -20°C. When reconstituted in sterile pH 7.4 0.01M PBS or diluent of antibody the antibody is stable for at least two weeks at 2-4 °C.
PubMed:	PubMed
Product Detail:	D-DIMER is A fibrin degradation fragment or product that is produced by the action of plasmin on fibrin in the clot dissolution process. SWISS: N/A Gene ID: N/A

Important Note:

This product as supplied is intended for research use only, not for use in human, therapeutic or diagnostic applications.


D-

二聚体是纤溶亢进的分子Marker之一。当凝血发生后，所形成的交联纤维在纤溶酶的作用下，两个D碎片交联成D-二聚体，它具有特殊的抗原性.作为DIC的重要早期诊断依据之一。

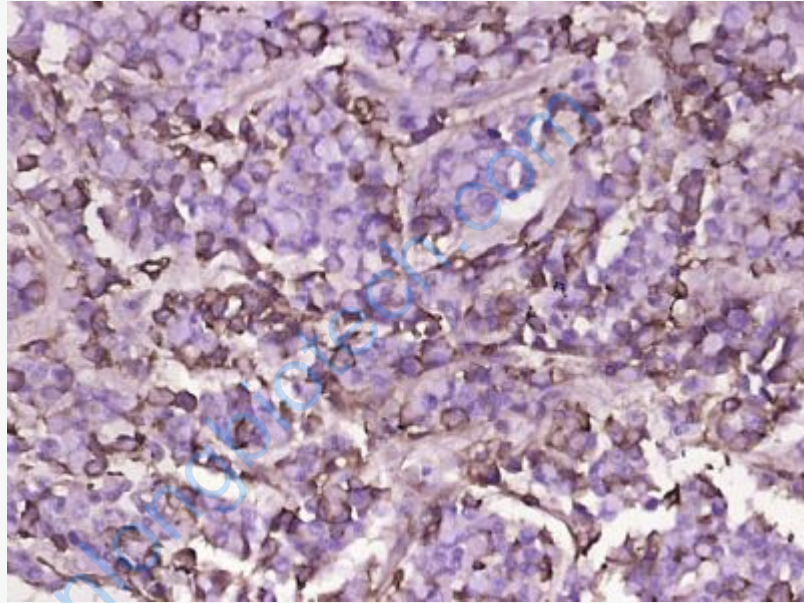

Picture:

Paraformaldehyde-fixed, paraffin embedded (human gastric carcinoma); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (D-DIMER) Polyclonal Antibody, Unconjugated (SL3514R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.

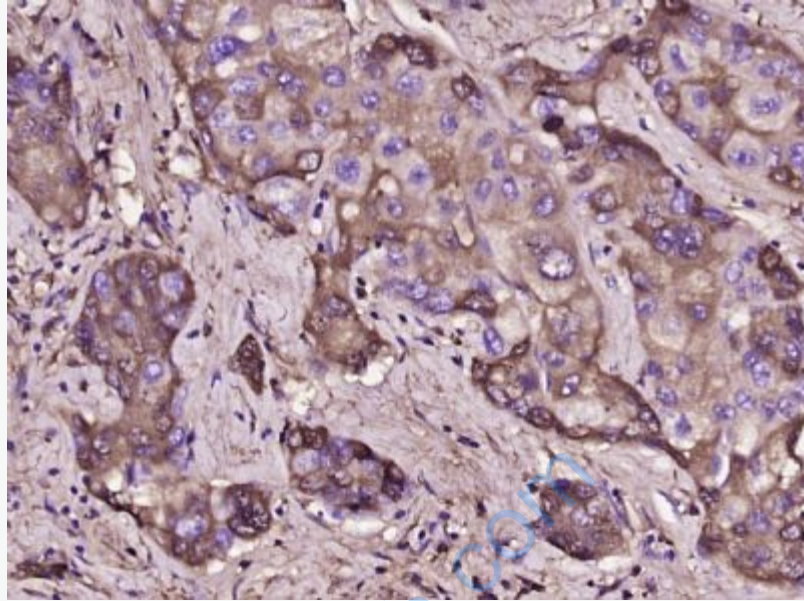

Paraformaldehyde-fixed, paraffin embedded (human liver carcinoma); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (D-DIMER) Polyclonal Antibody, Unconjugated (SL3514R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.