


Rabbit Anti-Phospholipase D2 antibody


SL8747R

Product Name:	Phospholipase D2
Chinese Name:	胆碱磷酸酯酶2抗体
Alias:	Choline phosphatase 2; EC 3.1.4.4; hPLD 2; hPLD2; mPLD2; Phosphatidylcholine hydrolyzing phospholipase D2; Phosphatidylcholine-hydrolyzing phospholipase D2; Phospholipase D2; PLD 2; PLD1C; Pld2; PLD2_HUMAN; Pldc; rPLD2.
Organism Species:	Rabbit
Clonality:	Polyclonal
React Species:	Human, Mouse, Rat, Dog, Cow, Horse, Rabbit, Sheep,
Applications:	ELISA=1:500-1000IHC-P=1:400-800IHC-F=1:400-800ICC=1:100-500IF=1:100-500 (Paraffin sections need antigen repair) not yet tested in other applications. optimal dilutions/concentrations should be determined by the end user.
Molecular weight:	106kDa
Cellular localization:	The cell membrane
Form:	Lyophilized or Liquid
Concentration:	lmg/ml
immunogen:	KLH conjugated synthetic peptide derived from human PLD2/Phospholipase D2:21-120/933
Lsotype:	IgG
Purification:	affinity purified by Protein A
Storage Buffer:	0.01M TBS(pH7.4) with 1% BSA, 0.03% Proclin300 and 50% Glycerol.
Storage:	Store at -20 °C for one year. Avoid repeated freeze/thaw cycles. The lyophilized antibody is stable at room temperature for at least one month and for greater than a year when kept at -20 °C. When reconstituted in sterile pH 7.4 0.01M PBS or diluent of antibody the antibody is stable for at least two weeks at 2-4 °C.
PubMed:	PubMed
Product Detail:	

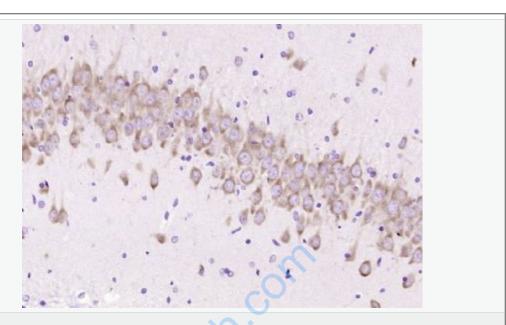

Paraformaldehyde-fixed, paraffin embedded (rat brain); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Phospholipase D2) Polyclonal Antibody, Unconjugated (SL8747R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.

Picture:

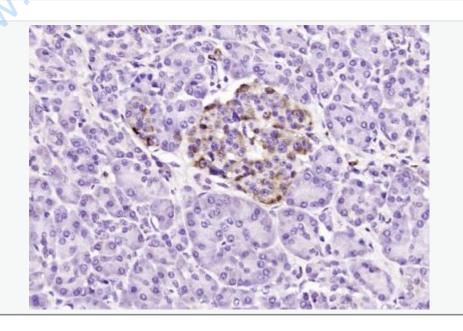

Paraformaldehyde-fixed, paraffin embedded (human pancreatic cancer); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Phospholipase D2) Polyclonal Antibody, Unconjugated (SL8747R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining. MANN SURPOROBIOTE CHI.CL