



SunLong Biotech Co.,LTD

Tel: 0086-571-56623320 Fax:0086-571-56623318

E-mail:sales@sunlongbiotech.com

www.sunlongbiotech.com

---

# **Human Regenerating Islet Derived Protein 3 Beta (REG3beta) ELISA Kit**

**96 Tests**

**Catalogue Number: SL2075Hu**

**Store all reagents at 2-8℃**

**Valied Period: six months**

**For samples:**

**Serum, plasma, cell culture supernatants, body fluid and tissue homogenate.**

**FOR RESEARCH USE ONLY !**

**NOT FOR THERAPEUTIC OR DIAGNOSTIC APPLICATIONS !**

**PLEASE READ THROUGH ENTIRE PROCEDURE BEFORE BEGINNING !**

# Human Regenerating Islet Derived Protein 3 Beta (REG3beta) ELISA Kit

**FOR RESEARCH USE ONLY**

---

---

## Drug Names

Generic Name: **Human Regenerating Islet Derived Protein 3 Beta (REG3beta) ELISA Kit**

## Purpose

Our Human Regenerating Islet Derived Protein 3 Beta (REG3beta) ELISA kit is to assay REG3beta levels in Human serum, plasma, culture media or any biological fluid.

## Principle

This ELISA kit uses Sandwich-ELISA as the method. The Microelisa stripplate provided in this kit has been pre-coated with an antibody specific to REG3beta. Standards or samples are added to the appropriate Microelisa stripplate wells and combined to the specific antibody. Then a Horseradish Peroxidase (HRP)-conjugated antibody specific for REG3beta is added to each Microelisa stripplate well and incubated. Free components are washed away. The TMB substrate solution is added to each well. Only those wells that contain REG3beta and HRP conjugated REG3beta antibody will appear blue in color and then turn yellow after the addition of the stop solution. The optical density (OD) is measured spectrophotometrically at a wavelength of 450 nm. The OD value is proportional to the concentration of REG3beta. You can calculate the concentration of REG3beta in the samples by comparing the OD of the samples to the standard curve.

## Materials provided with the kit

	Materials provided with the kit	96 determinations	Storage
1	User manual	1	R.T.
2	Closure plate membrane	2	R.T.
3	Sealed bags	1	R.T.
4	Microelisa stripplate	1	2-8°C

5	Standard: 1350pg/ml	0.5ml×1 bottle	2-8°C
6	Standard diluent	1.5ml×1 bottle	2-8°C
7	HRP-Conjugate reagent	6ml×1 bottle	2-8°C
8	Sample diluent	6ml×1 bottle	2-8°C
9	Chromogen Solution A	6ml×1 bottle	2-8°C
10	Chromogen Solution B	6ml×1 bottle	2-8°C
11	Stop Solution	6ml×1 bottle	2-8°C
12	wash solution	20ml (30X)×1bottle	2-8°C

## Sample preparation

### 1. Serum preparation

After collection of the whole blood, allow the blood to clot by leaving it undisturbed at room temperature. This usually takes 10-20 minutes. Remove the clot by centrifuging at 2,000-3,000 rpm for 20 minutes. If precipitates appear during reservation, the sample should be centrifugated again.

### 2. Plasma preparation

Collect the whole blood into tubes with anticoagulant (EDTA or citrate). After incubated at room temperature for 10-20 minutes, tubes are centrifugated for 20 min at 2,000-3,000 rpm. Collect the supernatant carefully as plasma samples. If precipitates appear during reservation, the sample should be centrifugated again.

### 3. Urine samples

Collect urine into aseptic tubes. Collect the supernatant carefully after centrifuging for 20 min at 2,000-3,000 rpm. If precipitates appear during reservation, the sample should be centrifugated again. The preparation procedure of cerebrospinal fluid and pleuroperitoneal fluid is the same as that of urine sample.

### 4. Cell samples

If you want to detect the secretions of cells, collect culture supernatant into aseptic tubes. Collect the supernatant carefully after centrifuging for 20 min at 2,000-3,000 rpm. If you want to detect intracellular components, dilute the cells to 1X100/ml with PBS (pH 7.2-7.4). The cells were destroyed to release intracellular components by repeated freezing and thawing. Collect the supernatant carefully after centrifuging for 20 min at 2,000-3,000 rpm. If precipitates appear during reservation, the sample should be centrifugated again.

## 5. Tissue samples

Tissue samples are cut, weighed, frozen in liquid nitrogen and stored at  $-80^{\circ}\text{C}$  for future use. The tissue samples were homogenized after adding PBS (pH 7.4). Samples should be operated at  $4^{\circ}\text{C}$ . Collect the supernatant carefully after centrifuging for 20 min at 2,000-3,000 rpm. Aliquot the supernatant for ELISA assay and future use.

### Notes:


1. Sample extraction and ELISA assay should be performed as soon as possible after sample collection. The samples should be extracted according to the relevant literature. If ELISA assay can not be performed immediately, samples can be stored at  $-20^{\circ}\text{C}$ . Repeated freeze-thaw cycles should be avoided.
2. Our kits can not be used for samples with  $\text{NaN}_3$  which can inhibit the activity of HRP.

## Procedure

### 1. Dilution of Standards

Dilute the standard by small tubes first, then pipette the volume of  $50\mu\text{l}$  from each tube to microplate well, each tube use two wells, total ten wells.

900pg/ml	Standard No.1	300 $\mu\text{l}$ Original Standard + 150 $\mu\text{l}$ Standard diluents
600pg/ml	Standard No.2	300 $\mu\text{l}$ Standard No.1 + 150 $\mu\text{l}$ Standard diluents
300pg/ml	Standard No.3	150 $\mu\text{l}$ Standard No.2 + 150 $\mu\text{l}$ Standard diluent
150pg/ml	Standard No.4	150 $\mu\text{l}$ Standard No.3 + 150 $\mu\text{l}$ Standard diluent
75pg/ml	Standard No.5	150 $\mu\text{l}$ Standard No.4 + 150 $\mu\text{l}$ Standard diluent

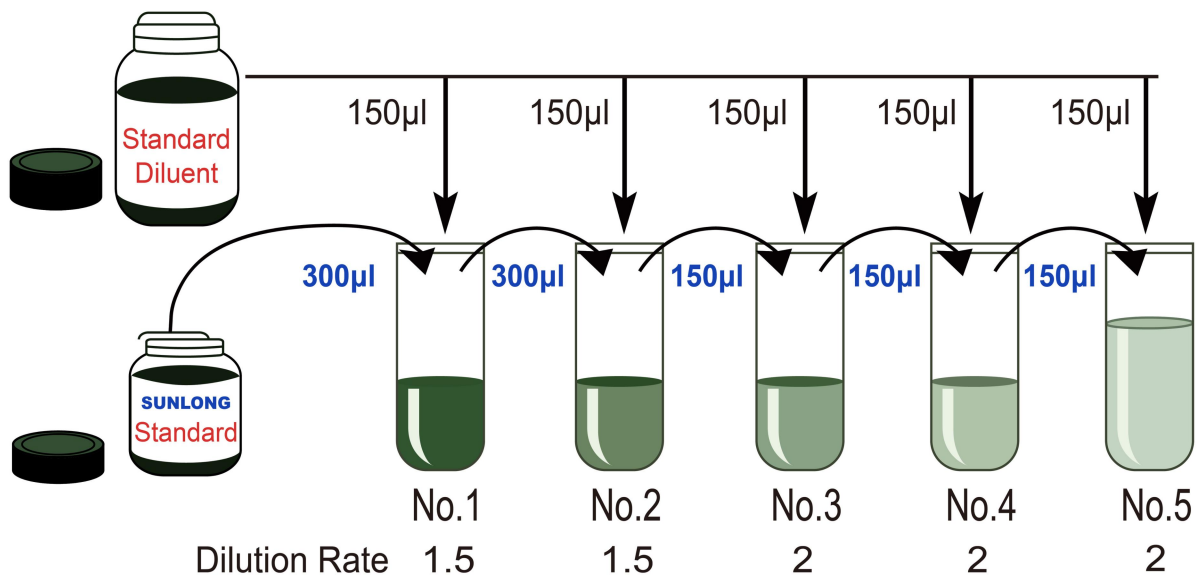

1350pg/ml	900pg/ml	600pg/ml	300pg/ml	150pg/ml	75pg/ml
-----------	----------	----------	----------	----------	---------

2. In the Microelisa stripplate, leave a well empty as blank control. In sample wells, 40µl Sample dilution buffer and 10µl sample are added (dilution factor is 5). Samples should be loaded onto the bottom without touching the well wall. Mix well with gentle shaking.
3. Incubation: incubate 30 min at 37°C after sealed with Closure plate membrane.
4. Dilution: dilute the concentrated washing buffer with distilled water (30 times for 96T and 20 times for 48T).
5. Washing: carefully peel off Closure plate membrane, aspirate and refill with the wash solution. Discard the wash solution after resting for 30 seconds. Repeat the washing procedure for 5 times.
6. Add 50 µl HRP-Conjugate reagent to each well except the blank control well.
7. Incubation as described in Step 3.
8. Washing as described in Step 5.
9. Coloring: Add 50 µl Chromogen Solution A and 50 µl Chromogen Solution B to each well, mix with gently shaking and incubate at 37°C for 15 minutes. Please avoid light during coloring.
10. Termination: add 50 µl stop solution to each well to terminate the reaction. The color in the well should change from blue to yellow.
11. Read absorbance O.D. at 450nm using a Microtiter Plate Reader. The OD value of the blank control well is set as zero. Assay should be carried out within 15 minutes after adding stop solution.


### **Notes:**

1. Store the kit at 4° C upon receipt. The kit should be equilibrated to room temperature before the assay. Remove any unneeded strips from Human REG3beta Antibody-Coated plate, reseal them in zip-lock foil and keep at 4° C.
2. Precipitates may appear in concentrated washing buffer. Please heat the buffer to dissolve all the precipitates, which will not affect the results.
3. Accurate pipette should be used to avoid experimental error. Samples should be added to the Microplate in less than 5 minutes. If a large number of samples are included, multiple channel pipette is recommended.

4. Standard curve should be included in every assay. Replicate wells are recommended. If the OD value of the sample is greater than the first well of standards, please dilute the sample (n times) before test. When calculating the original REG3beta concentration, please multiply the total dilution factor (XnX5).
5. In order to avoid cross-contamination, Microplate sealers are for one-time use only.
6. Please keep Substrate away from light.
7. All the operation should be accordance with the manufacturer's instructions strictly. The results determined by the Microtiter Plate Reader.
8. All the samples, washing buffer and wastes should be treated as infectious agents.
9. Reagents from different lots should not be mixed.

## Calculation of Results

Known concentrations of Human REG3beta Standard and its corresponding reading OD is plotted on the log scale (x-axis) and the log scale (y-axis) respectively. The concentration of Human REG3beta in sample is determined by plotting the sample's O.D. on the Y-axis. The original concentration is calculated by multiplying the dilution factor.

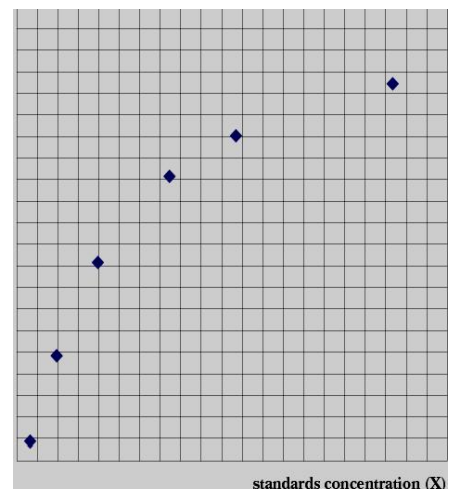

This diagram is for reference only

## Precision

Intra-assay Precision (Precision within an assay): 3 samples with low, middle and high level Human REG3beta were tested 20 times on one plate, respectively.

Inter-assay Precision (Precision between assays): 3 samples with low, middle and high level Human REG3beta were tested on 3 different plates, 8 replicates in each plate.

$$CV(\%) = SD/\text{mean} \times 100$$

Intra-Assay: CV<10%

Inter-Assay: CV<12%

### **Assay range**

16pg/ml-1000pg/ml

### **Sensitivity**

6.5pg/ml

### **Storage and validity**

1. Storage: 2-8°C.
2. Duration: 6 months