



SunLong Biotech Co.,LTD

Tel: 0086-571-56623320 Fax:0086-571-56623318

E-mail:sales@sunlongbiotech.com

www.sunlongbiotech.com

Human Oxidized Low Density Lipoprotein antibody IgG (OxLDL-IgG) ELISA Kit

96 Tests

Catalogue Number:SL3570Hu

Store all reagents at 2-8℃

Validity Period: six months

For samples:

In serum, plasma, culture media or any biological fluid.

FOR RESEARCH USE ONLY !

NOT FOR THERAPEUTIC OR DIAGNOSTIC APPLICATIONS !

PLEASE READ THROUGH ENTIRE PROCEDURE BEFORE BEGINNING !

Human Oxidized Low Density Lipoprotein antibody IgG (OxLDL-IgG) ELISA Kit

FOR RESEARCH USE ONLY

Drug Names

Generic Name: **Human Oxidized Low Density Lipoprotein antibody IgG (OxLDL-IgG) ELISA Kit**

Purpose

Our Human Oxidized Low Density Lipoprotein antibody IgG (OxLDL-IgG) ELISA kit is to assay OxLDL-IgG levels in Human serum, plasma, culture media or any biological fluid.

Principle

This ELISA kit uses Sandwich-ELISA as the method. The Microelisa stripplate provided in this kit has been pre-coated with an antigen specific to OxLDL-IgG. Standards or samples are added to the appropriate Microelisa stripplate wells and combined to the specific antigen. Then a Horseradish Peroxidase (HRP)-conjugated antigen specific for OxLDL-IgG is added to each Microelisa stripplate well and incubated. Free components are washed away. The TMB substrate solution is added to each well. Only those wells that contain OxLDL-IgG and HRP conjugated OxLDL antigen will appear blue in color and then turn yellow after the addition of the stop solution. The optical density (OD) is measured spectrophotometrically at a wavelength of 450 nm. The OD value is proportional to the concentration of OxLDL-IgG. You can calculate the concentration of OxLDL-IgG in the samples by comparing the OD of the samples to the standard curve.

Materials provided with the kit

	Materials provided with the kit	96 determinations	Storage
1	User manual	1	R.T.
2	Closure plate membrane	2	R.T.
3	Sealed bags	1	R.T.
4	Microelisa stripplate	1	2-8°C
5	Standard: 360 μ g/L	0.5ml×1 bottle	2-8°C
6	Standard diluent	1.5ml×1 bottle	2-8°C
7	HRP-Conjugate reagent	6ml×1 bottle	2-8°C

8	Sample diluent	6ml×1 bottle	2-8°C
9	Chromogen Solution A	6ml×1 bottle	2-8°C
10	Chromogen Solution B	6ml×1 bottle	2-8°C
11	Stop Solution	6ml×1 bottle	2-8°C
12	wash solution	20ml (30X)×1bottle	2-8°C

Sample preparation

1. Serum preparation

After collection of the whole blood, allow the blood to clot by leaving it undisturbed at room temperature. This usually takes 10-20 minutes. Remove the clot by centrifuging at 2,000-3,000 rpm for 20 minutes. If precipitates appear during reservation, the sample should be centrifugated again.

2. Plasma preparation

Collect the whole blood into tubes with anticoagulant (EDTA or citrate). After incubated at room temperature for 10-20 minutes, tubes are centrifugated for 20 min at 2,000-3,000 rpm. Collect the supernatant carefully as plasma samples. If precipitates appear during reservation, the sample should be centrifugated again.

3. Urine samples

Collect urine into aseptic tubes. Collect the supernatant carefully after centrifuging for 20 min at 2,000-3,000 rpm. If precipitates appear during reservation, the sample should be centrifugated again. The preparation procedure of cerebrospinal fluid and pleuroperitoneal fluid is the same as that of urine sample.

4. Cell samples

If you want to detect the secretions of cells, collect culture supernatant into aseptic tubes. Collect the supernatant carefully after centrifuging for 20 min at 2,000-3,000 rpm. If you want to detect intracellular components, dilute the cells to 1×10^6 /ml with PBS (pH 7.2-7.4). The cells were destroyed to release intracellular components by repeated freezing and thawing. Collect the supernatant carefully after centrifuging for 20 min at 2,000-3,000 rpm. If precipitates appear during reservation, the sample should be centrifugated again.

5. Tissue samples

Tissue samples are cut, weighed, frozen in liquid nitrogen and stored at -80°C for future

use. The tissue samples were homogenized after adding PBS (pH 7.4). Samples should be operated at 4°C. Collect the supernatant carefully after centrifuging for 20 min at 2,000-3,000 rpm. Aliquot the supernatant for ELISA assay and future use.

Notes:

1. Sample extraction and ELISA assay should be performed as soon as possible after sample collection. The samples should be extracted according to the relevant literature. If ELISA assay can not be performed immediately, samples can be stored at -20 °C .Repeated freeze-thaw cycles should be avoided.
2. Our kits can not be used for samples with NaN₃ which can inhibit the activity of HRP.


Assay procedure

Procedure

1. Dilution of Standards

Dilute the standard by small tubes first, then pipette the volume of 50ul from each tube to microplate well, each tube use two wells ,total ten wells.

240 μ g/L	Standard No.1	300μl Original Standard + 150μl Standard diluents
160 μ g/L	Standard No.2	300μl Standard No.1 + 150μl Standard diluents
80 μ g/L	Standard No.3	150μl Standard No.2 + 150μl Standard diluent
40 μ g/L	Standard No.4	150μl Standard No.3 + 150μl Standard diluent
20 μ g/L	Standard No.5	150μl Standard No.4 + 150μl Standard diluent

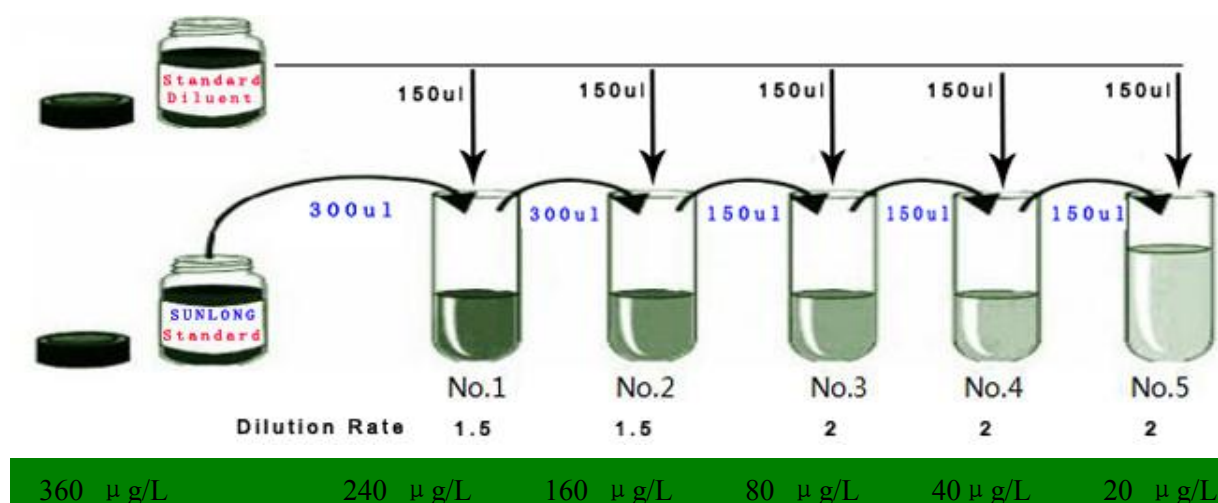

2. In the Microelisa stripplate, leave a well empty as blank control. In sample wells, 40μl Sample dilution buffer and 10μl sample are added (dilution factor is 5). Samples should be

loaded onto the bottom without touching the well wall. Mix well with gentle shaking.

3. Incubation: incubate 30 min at 37°C after sealed with Closure plate membrane.
4. Dilution: dilute the concentrated washing buffer with distilled water (30 times for 96T and 20 times for 48T).
5. Washing: carefully peel off Closure plate membrane, aspirate and refill with the wash solution. Discard the wash solution after resting for 30 seconds. Repeat the washing procedure for 5 times.
6. Add 50 µl HRP-Conjugate reagent to each well except the blank control well.
7. Incubation as described in Step 3.
8. Washing as described in Step 5.
9. Coloring: Add 50 µl Chromogen Solution A and 50 µl Chromogen Solution B to each well, mix with gently shaking and incubate at 37°C for 15 minutes. Please avoid light during coloring.
10. Termination: add 50 µl stop solution to each well to terminate the reaction. The color in the well should change from blue to yellow.
11. Read absorbance O.D. at 450nm using a Microtiter Plate Reader. The OD value of the blank control well is set as zero. Assay should be carried out within 15 minutes after adding stop solution.

Notes:


1. Store the kit at 4 ° C upon receipt. The kit should be equilibrated to room temperature before the assay. Remove any unneeded strips from Human OxLDL antigen-Coated plate, reseal them in zip-lock foil and keep at 4 ° C.
2. Precipitates may appear in concentrated washing buffer. Please heat the buffer to dissolve all the precipitates, which will not affect the results.
3. Accurate pipette should be used to avoid experimental error. Samples should be added to the Microplate in less than 5 minutes. If a large number of samples are included, multiple channel pipette is recommended.
4. Standard curve should be included in every assay. Replicate wells are recommended. If the OD value of the sample is greater than the first well of standards, please dilute the

sample (n times) before test. When calculating the original OxLDL-IgG concentration, please multiply the total dilution factor (XnX5).

5. In order to avoid cross-contamination, Microplate sealers are for one-time use only.
6. Please keep Substrate away from light.
7. All the operation should be accordance with the manufacturer's instructions strictly. The results determined by the Microtiter Plate Reader.
8. All the samples, washing buffer and wastes should be treated as infectious agents.
9. Reagents from different lots should not be mixed.

Calculation of Results

Known concentrations of Human OxLDL-IgG Standard and its corresponding reading OD is plotted on the log scale (x-axis) and the log scale (y-axis) respectively. The concentration of Human OxLDL-IgG in sample is determined by plotting the sample's O.D. on the Y-axis. The original concentration is calculated by multiplying the dilution

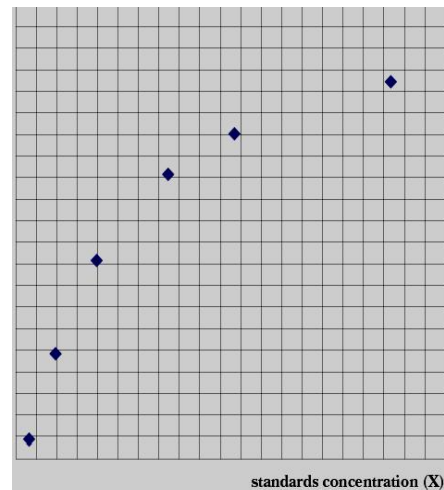

This diagram is for reference only

Precision

Intra-assay Precision (Precision within an assay): 3 samples with low, middle and high level Human OxLDL-IgG were tested 20 times on one plate, respectively.

Inter-assay Precision (Precision between assays): 3 samples with low, middle and high level Human OxLDL-IgG were tested on 3 different plates, 8 replicates in each plate.

$$CV(\%) = SD/\text{mean} \times 100$$

Intra-Assay: $CV < 10\%$

Inter-Assay: $CV < 12\%$

Assay range

5 μ g/L - 330 μ g/L

Sensitivity:

1 μ g/L

Storage and validity

1. Storage: 2-8°C.
2. Duration: 6 months