# Active Glial Cell Line Derived Neurotrophic Factor (GDNF) Instruction Manual

## SBPA026Hu01

### Homo sapiens (Human)

**Buffer Formulation**20mM Tris, 150mM NaCl, pH8.0, containing 1mM EDTA, 1mM DTT, 0.01% SKL, 5% Trehalose and Proclin300.

**Traits** Freeze-dried powder

Purity > 95% Isoelectric Point 9.4

**Applications** Cell culture; Activity Assays.

#### **ACTIVITY TEST**

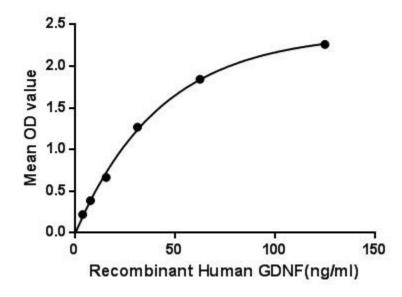


Figure 1. The binding activity of GDNF with GFRa2.

#### **USAGE**

Reconstitute in 20mM Tris, 150mM NaCl (PH8.0) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

#### **STORAGE**

Avoid repeated freeze/thaw cycles. Store at 2-8°C for one month. Aliquot and store at -80°C for 12 months.

#### **STABILITY**

The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

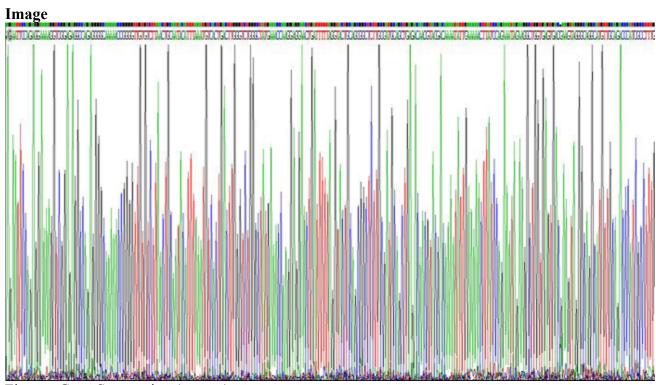


Figure . Gene Sequencing (extract)

## **Image**

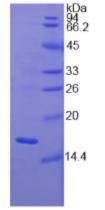


Figure. SDS-PAGE

# [IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.