Active Matrix Metalloproteinase 2 (MMP2) Instruction Manual

SBPA067Hu02

Homo sapiens (Human)

Buffer Formulation PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5%


Trehalose and Proclin300.

Traits Freeze-dried powder

Purity > 97% Isoelectric Point 5.2

Applications Cell culture; Activity Assays.

ACTIVITY TEST

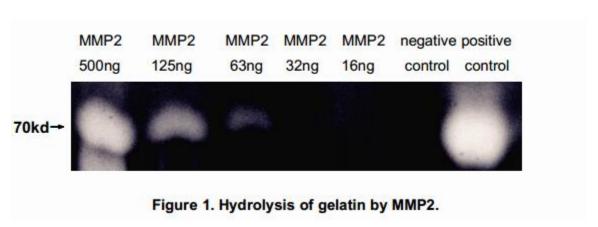

Mechanism: MMP2 is a zinc-dependent enzymes capable of cleaving components of the extracellular matrix, which belongs to the matrix metalloproteinase (MMP) family. It is a gelatinase A, 72kDa type IV collagenase which can hydrolyze gelatin under certain conditions. Gelatin zymography is mainly used for the detection of the gelatinases, MMP-2 and MMP-9 and It is extremely sensitive because levels of 10pg of MMP-2 can already be detected. Briefly, various concentrations of MMP2 (500ng, 125ng, 63ng, 32ng, 16ng) were denatured by SDS loading buffer, electrophoresed through sodium dodecylsulphate-polyacrylamide gel (SDS-PAGE; 10% gels) containing gelatin (1mg/mL) with nonreducing conditions. After renaturation, incubation and CCB-stained, active MMP2 would hydrolyze gelatin nearby, which was indicated by the white binds on the gel. In this experiment we use heat-denatured MMP2 protein as negative control, and blood sample as positive control.

Result 1: Gelatin hydrolysis by recombinant human MMP2 was shown in figure 1.

USAGE


Reconstitute in PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

STORAGE

Avoid repeated freeze/thaw cycles. Store at 2-8°C for one month. Aliquot and store at -80°C for 12 months.

STABILITY

The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

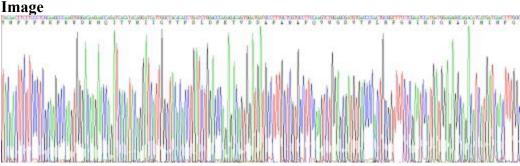

SDS-PAGE Image

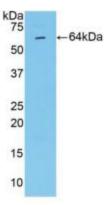


Figure. Western Blot; Sample: Recombinant MMP2, Human.

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.