Active Janus Kinase 3 (JAK3) Instruction Manual

SBPF358Hu01

Homo sapiens (Human)

Buffer Formulation 20mM Tris, 150mM NaCl, pH8.0, containing 1mM EDTA,

1mM DTT, 0.01% SKL, 5% Trehalose and Proclin300.

Traits Freeze-dried powder

Purity > 97% Isoelectric Point 8.5

Applications Cell culture; Activity Assays.

ACTIVITY TEST

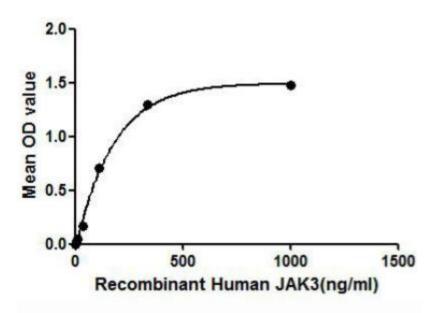


Figure 1. The binding activity of JAK3 with STAM2.

JAK3 (Tyrosine-protein kinase JAK3) is a non-receptor tyrosine kinase that involved in various processes such as cell growth, development, or differentiation. Besides, STAM2 has been identified as an interactor of JAK3, thus a binding ELISA assay was conducted to detect the interaction of recombinant human JAK3 and recombinant human STAM2. Briefly, JAK3 were diluted serially in PBS, with 0.01%BSA (pH 7.4). Duplicate samples of 100uL were then transferred to STAM2-coated microtiter wells and incubated for 2h at 37°C. Wells were washed with PBST and incubated for 1h with anti-JAK3 pAb, then aspirated and washed 3 times. After incubation with HRP labelled secondary antibody, wells were aspirated and washed 3 times. With the addition of substrate solution, wells

were incubated 15-25 minutes at 37°C. Finally, add 50µL stop solution to the wells and read at 450nm immediately. The binding activity of of JAK3 and STAM2 was shown in Figure 1, and this effect was in a dose dependent manner.

USAGE


Reconstitute in 20mM Tris, 150mM NaCl (pH8.0) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

STORAGE

Avoid repeated freeze/thaw cycles. Store at 2-8°C for one month. Aliquot and store at -80°C for 12 months.

STABILITY

The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

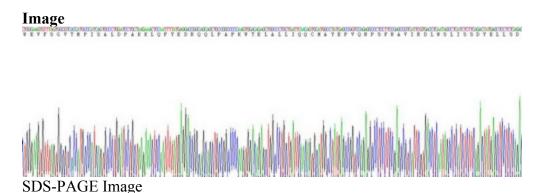

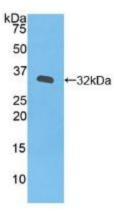


Figure. Western Blot; Sample: Recombinant JAK3, Human.

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.