

# Recombinant Matrix Metalloproteinase 2 (MMP2) Instruction Manual

**SIPA067Hu02**

**Homo sapiens (Human)**

<b>Source</b>	Prokaryotic expression
<b>Host</b>	E.coli
<b>Endotoxin Level</b>	<1.0EU per 1µg (determined by the LAL method)
<b>Subcellular Location</b>	Membrane, Nucleus, Extracellular matrix
<b>Predicted Molecular Mass</b>	63.3kDa
<b>Accurate Molecular Mass</b>	63kDa(Analysis of differences refer to the manual)
<b>Residues &amp; Tags</b>	Tyr110~Cys660 with N-terminal His Tag
<b>Buffer Formulation</b>	20mM Tris, 150mM NaCl, pH8.0, containing 0.01% SKL, 5% Trehalose.
<b>Traits</b>	Freeze-dried powder
<b>Purity</b>	> 97%
<b>Isoelectric Point</b>	5.2
<b>Applications</b>	Positive Control; Immunogen; SDS-PAGE; WB.

## SEQUENCE

```
Y NFFPRKPKWD KNQITYRIIG YTPDLDPETV DDAFARAFQV  
WSDVTPLRFS RIHDGEADIM INFGRWEHGD GYPFDGKDGL LAHAFAPGTG  
VGGDSHFDDD ELWTLGEGQV VRVKYGNADG EYCKFPFLFN GKEYNSCTDT  
GRSDGFLWCS TTYNFEKDGK YGFCPHEALF TMGGNAEQQP CKFPFRFQGT  
SYDSCTTEGR TDGYRWCGTT EDYDRDKKYG FCPETAMSTV GGNSEGAPCV  
FPFTFLGNKY ESCTSAGRSD GKMWCATTAN YDDRKKGFC PDQGYSLFLV  
AAHEFGHAMG LEHSQDPGAL MAPIYTYTKN FRLSQDDIKG IQELYGASPD  
IDLGTGPTPT LGPVTPEICK QDIVFDGIAQ IRGEIFFFKD RFIWRTVTPR  
DKPMGPLLVA TFWPELPEKI DAVYEAPQEE KAVFFAGNEY WIYSASTLER  
GYPKPLTSLG LPPDVQRVDA AFNWSKNKKT YIFAGDKFWR YNEVKKKMDP  
GFPKLIADAW NAIPDNLDAV VDLQGGGHSY FFKGAYYLKL ENQSLKSVKF  
GSIKSDWLGC
```

## USAGE

Reconstitute in 20mM Tris, 150mM NaCl (pH8.0) to a concentration of 0.1-1.0 mg/mL. Do not vortex.


## STORAGE

Avoid repeated freeze/thaw cycles. Store at 2-8°C for one month. Aliquot and store at -80°C for 12 months.

## STABILITY

The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## Image

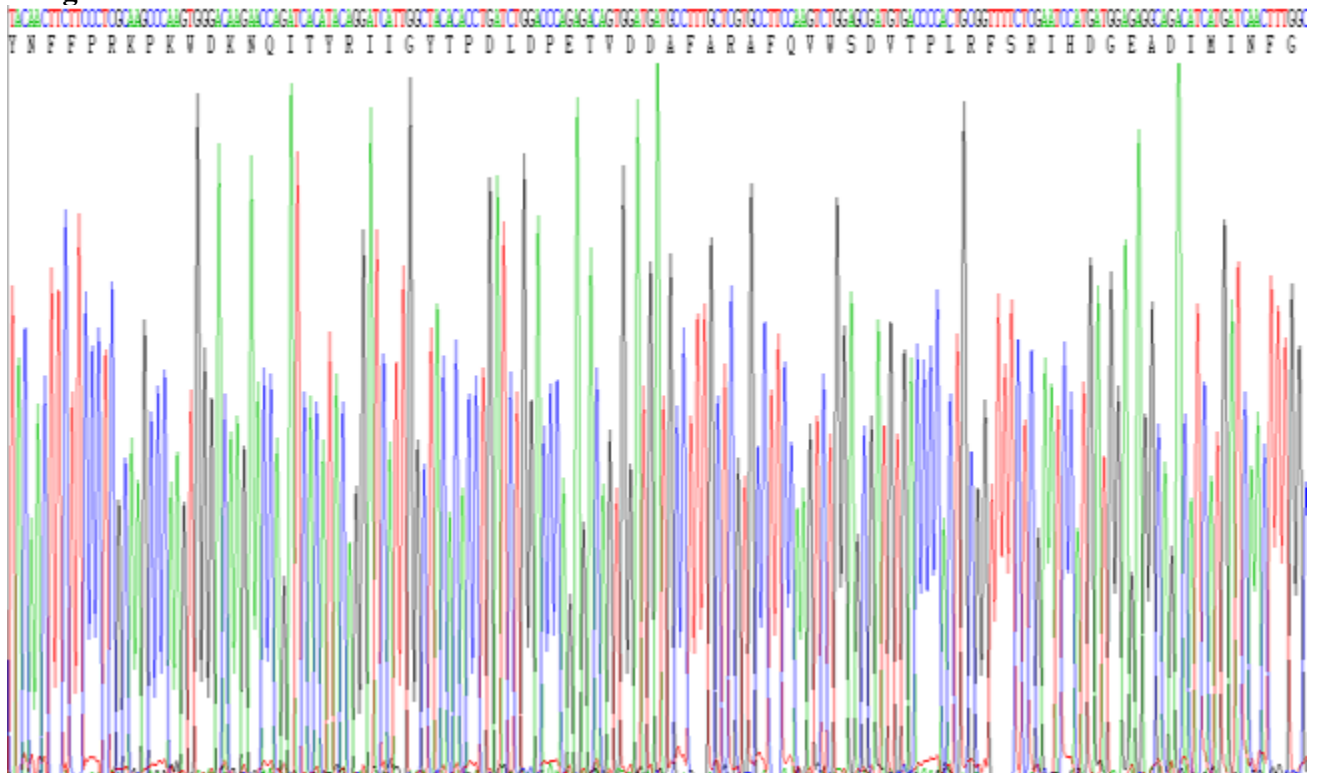

SDS-PAGE Image

**[IMPORTANT NOTE]**

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.