

Recombinant Centromere Protein I (CENPI) Instruction Manual

SIPA337Hu01

Homo sapiens (Human)

Source	Prokaryotic expression
Host	E.coli
Endotoxin Level	<1.0EU per 1µg (determined by the LAL method)
Subcellular Location	Nucleus, Chromosome
Predicted Molecular Mass	90.4kDa
Accurate Molecular Mass	90kDa(Analysis of differences refer to the manual)
Residues & Tags	Met1~Tyr756 with N-terminal His Tag
Buffer Formulation	PBS, pH7.4, containing 0.01% SKL, 5% Trehalose.
Traits	Freeze-dried powder
Purity	> 95%
Isoelectric Point	9.0
Applications	Positive Control; Immunogen; SDS-PAGE; WB.

SEQUENCE

```
MSPQKRVKNVQAQNRTSQGSSTFQTTL SAKWKVKQDPSNSKNISKHGQNNPVGDIYEHADDQAEEDALQMAVGYFEKGP IKASQNKDKTL  
EKHLKTVENVAWKNGLASEE IDILLNIALSGKFGNAVNTRILKCMIPATVISEDVVKAVSWLCVKGKCSGSTKVL FYRWLVAMFDFID  
RKEQINLLYGFFASLQDDALCPYVCHLLYLL TKKENVKPFRVRKLLDLQAKMGMQPHLQALL SLYKFFAPALISVSLPVRKKIYFKN  
SENLWKTALLAVKQRNRGSPPEPLKMLGPANVRPLKRWNSLSVIPVLNSSSYTKECGKKEMSLSDCLNRSGSFPLEQLQSFQQLLQ  
NIHCLELPSQMGSVLNNSLLLHYINCVRDEPVLLRFYYWLSQTLQEECIWKVNNYEHGKFTNF LDTIIRAECFLQEGFYSCFAFLY  
KSLPLWDGLCCRSQFLQLVSWIPFSSSFSEVKPLLFDHLAQLFFTSTIYFKCSVLQSLKELLQNWLLWLSMDIHMKPVTNSPLETTLGG  
SMNSVSKLIHYVGLSTTAMRLESNNTFLLHF ILDFYEKVCDIYINYNLPLVLFPPGIFYSALLSLDTSILNQLCFIMHRYRKNLTA  
AKKNELVQKTKSEFNFSKTYQEFNHYL TSMVGLWTSKPFKGKIYIDPEILEKTGVAEYKNSLNVVHHP SFLSYAVSFLQESPEER  
TVNVSSIRGKKWSWYLDYLF SQGLQGLKLFIRSSVHHSIPRAEGINCNNQY
```

USAGE

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

STORAGE

Avoid repeated freeze/thaw cycles. Store at 2-8°C for one month. Aliquot and store at -80°C for 12 months.

STABILITY

The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

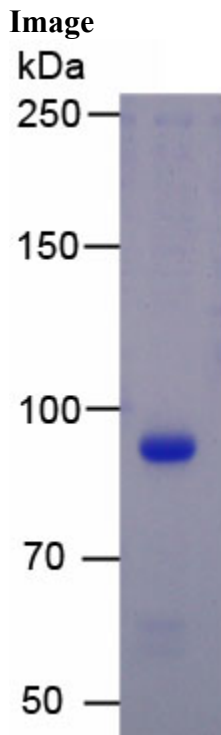


Figure. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.