

# Recombinant Collagen Type I Alpha 1 (COL1a1) Instruction Manual

## SIPA345Mu01

### Mus musculus (Mouse)

<b>Source</b>	Prokaryotic expression
<b>Host</b>	E.coli
<b>Endotoxin Level</b>	<1.0EU per 1µg (determined by the LAL method)
<b>Subcellular Location</b>	Secreted
<b>Predicted Molecular Mass</b>	26.8kDa
<b>Accurate Molecular Mass</b>	33kDa(Analysis of differences refer to the manual)
<b>Residues &amp; Tags</b>	Lys1225~Val1453 with N-terminal His Tag
<b>Buffer Formulation</b>	PBS, pH7.4, containing 0.05% SKL, 1mM DTT, 5% Trehalose and Proclin300.
<b>Traits</b>	Freeze-dried powder
<b>Purity</b>	> 95%
<b>Isoelectric Point</b>	5.7
<b>Applications</b>	Positive Control; Immunogen; SDS-PAGE; WB.

### SEQUENCE

```
                KSLSQQ IENIRSPEGS RKNPARTCRD
LKMCHSDWKS GEYWIDPNQG CNLDAIKVYC NMETGQTCVF PTQPSVPQKN
WYISPMPKEK KHWVFGESMT DGFPFEYGSE GSDPADVAIQ LTFLRLMSTE
ASQNITYHCK NSVAYMDQQT GNLKKALLLQ GSNEIELRGE GNSRFTYSTL
VDGCTSHTGT WGKTVIEYKT TKTSRLPIID VAPLDIGAPD QEFGLDIGPA
CFV
```

```
                KSLSQQ IENIRSPEGS RKNPARTCRD
LKMCHSDWKS GEYWIDPNQG CNLDAIKVYC NMETGQTCVF PTQPSVPQKN
WYISPMPKEK KHWVFGESMT DGFPFEYGSE GSDPADVAIQ LTFLRLMSTE
ASQNITYHCK NSVAYMDQQT GNLKKALLLQ GSNEIELRGE GNSRFTYSTL
VDGCTSHTGT WGKTVIEYKT TKTSRLPIID VAPLDIGAPD QEFGLDIGPA
CFV
```

### USAGE

Reconstitute in 10mM PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

## **STORAGE**

Avoid repeated freeze/thaw cycles. Store at 2-8°C for one month. Aliquot and store at -80°C for 12 months.

## **STABILITY**

The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## **Image**

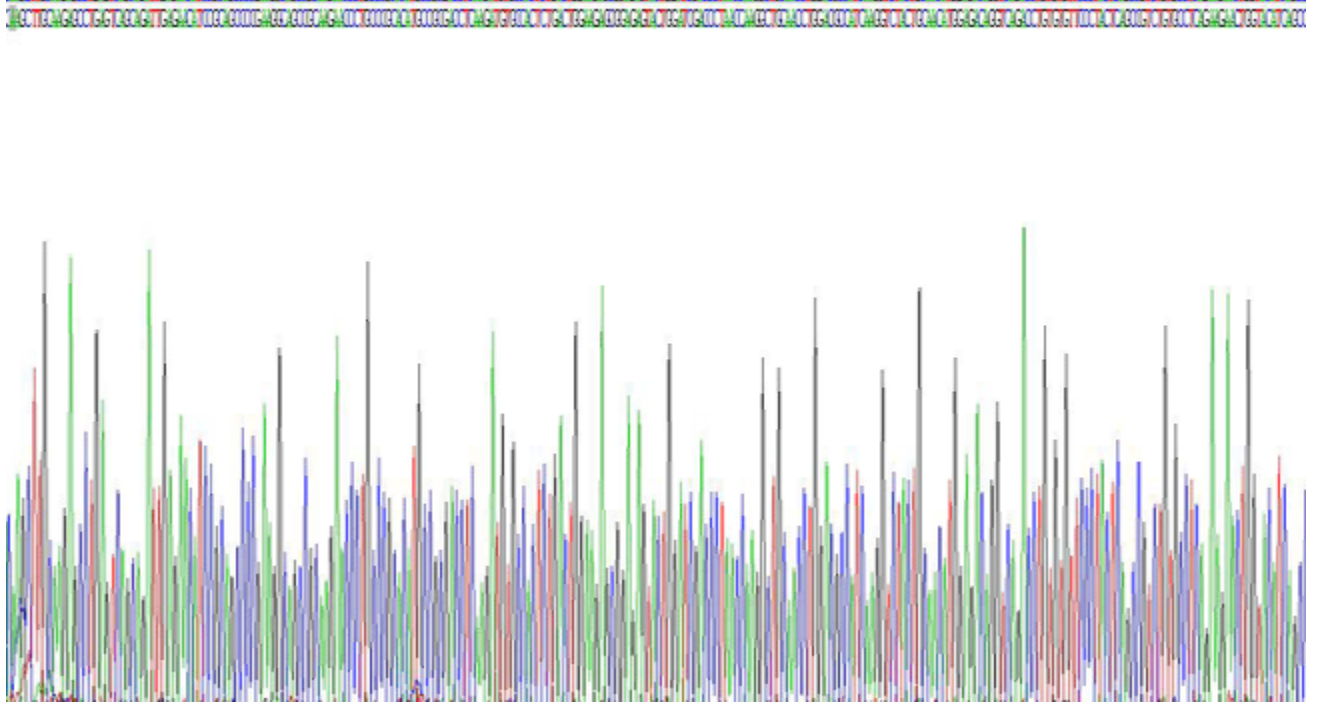


Figure. SDS-PAGE

## **[IMPORTANT NOTE]**

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.