# Recombinant Heat Shock 70kDa Protein 1A (HSPA1A) Instruction Manual

# **SIPB053B001**

**Bos taurus; Bovine (Cattle)** 

**Source** Prokaryotic expression

**Host** E.coli

Endotoxin Level <1.0EU per 1µg (determined by the LAL method)

Subcellular LocationSecretedPredicted Molecular Mass71.5kDa

Accurate Molecular Mass 70kDa(Analysis of differences refer to the manual)

Residues & Tags Met1~Asp641 with N-terminal His Tag

**Buffer Formulation** PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5%

Trehalose and Proclin300.

**Traits** Freeze-dried powder

Purity > 90% Isoelectric Point 5.9

**Applications** Positive Control; Immunogen; SDS-PAGE; WB.

## **SEQUENCE**

MAKNMAIGID LGTTYSCVGV FQHGKVEIIA NDQGNRTTPS YVAFTDTERL IGDAAKNQVA LNPQNTVFDA KRLIGRKFGD PVVQSDMKHW PFRVINDGDK PKVQVSYKGE TKAFYPEEIS SMVLTKMKEI AEAYLGHPVT NAVITVPAYF NDSQRQATKD AGVIAGLNVL RIINEPTAAA IAYGLDRTGK GERNVLIFDL GGGTFDVSIL TIDDGIFEVK ATAGDTHLGG EDFDNRLVNH FVEEFKRKHK KDISQNKRAV RRLRTACERA KRTLSSSTQA SLEIDSLFEG IDFYTSITRA RFEELCSDLF RSTLEPVEKA LRDAKLDKAQ IHDLVLVGGS TRIPKVQKLL QDFFNGRDLN KSINPDEAVA YGAAVQAAIL MGDKSENVQD LLLLDVAPLS LGLETAGGVM TALIKRNSTI PTKQTQIFTT YSDNQPGVLI QVYEGERAMT RDNNLLGRFE LSGIPPAPRG VPQIEVTFDI DANGILNVTA TDKSTGKANK ITITNDKGRL SKEEIERMVQ EAEKYKAEDE VQRERVSAKN ALESYAFNMK SAVEDEGLKG KISEADKKKV LDKCQEVISW LDANTLAEKD EFEHKRKELE OVCNPIISRL YOGAGGPGAG GFGAOGPKGG SGSGPTIEEV D

#### USAGE

Reconstitute in 10mM PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex.


#### **STORAGE**

Avoid repeated freeze/thaw cycles. Store at 2-8°C for one month. Aliquot and store at -80°C for 12 months.

## **STABILITY**

The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

#### **Image**

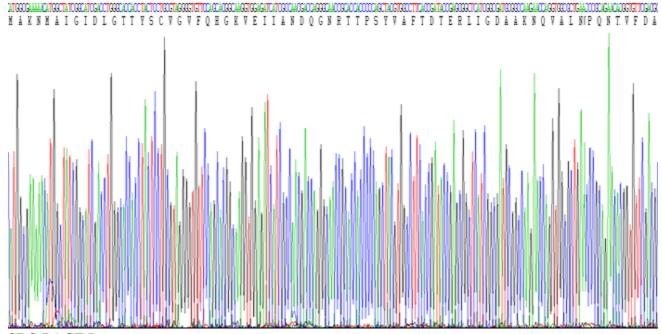

SDS-PAGE Image

# [IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.