

# Recombinant Polypyrimidine Tract Binding Protein 1 (PTBP1)

## Instruction Manual

**SIPC458Hu01**

**Homo sapiens (Human)**

<b>Source</b>	Prokaryotic expression
<b>Host</b>	E.coli
<b>Endotoxin Level</b>	<1.0EU per 1µg (determined by the LAL method)
<b>Subcellular Location</b>	Nucleus
<b>Predicted Molecular Mass</b>	60.9kDa
<b>Accurate Molecular Mass</b>	61kDa(Analysis of differences refer to the manual)
<b>Residues &amp; Tags</b>	Met1~Ile531 with N-terminal His Tag
<b>Buffer Formulation</b>	100mMNaHCO <sub>3</sub> , 500mMNaCl, pH8.3, containing 1mM EDTA, 1mM DTT, 0.01% SKL, 5% Trehalose and Proclin300.
<b>Traits</b>	Freeze-dried powder
<b>Purity</b>	> 90%
<b>Isoelectric Point</b>	9.2
<b>Applications</b>	Positive Control; Immunogen; SDS-PAGE; WB.

### SEQUENCE

MDGIVPDI AV GTKR GSDE LF STCV NGPFI MSSNSASAAN GNDSKKFKGD  
SRSAGVPSRV IHIRKLPIDV TEGEVISLGL PFGKVTNLLM LKGKNQAFIE  
MNTEEAANTM VNYYTSVTPV LRQPIYIQF SNHKELKTD S SPNQARAQAA  
LQAVNSVQSG NLALAASAAA VDAGMAMAGQ SPVLRIIVEN LFYPVTLDVL  
HQIFS KFGTV LKIITFTKNN QFQALLQYAD PVSAQHAKLS LDGQNIYNAC  
CTLRIDFSKL TSLNVKYNND KSRDYTRPD L PSGDSQPSLD QTMAAAFGLS  
VPNVHGALAP LAIPSAAAAA AAAGRIAIPG LAGAGNSVLL VSNLNPERVT  
PQSLFILFGV YGDVQRVKIL FNKKENALVQ MADGNQAQLA MSHLNGHKLH  
GKPIRITLSK HQNVQLPREG QEDQGLTKDY GNSPLHRFKK PGSKNFQNIF  
PPSATLHLSN IPPSVSEEDL KVLFSSNGGV VKGFKFQKD RKMALIQMGS  
VEEAVQALID LHNHDLGENH HLRVSFSKST I

### USAGE

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-1.0 mg/mL. Do not vortex.


## STORAGE

Avoid repeated freeze/thaw cycles. Store at 2-8°C for one month. Aliquot and store at -80°C for 12 months.

## STABILITY

The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

### Image

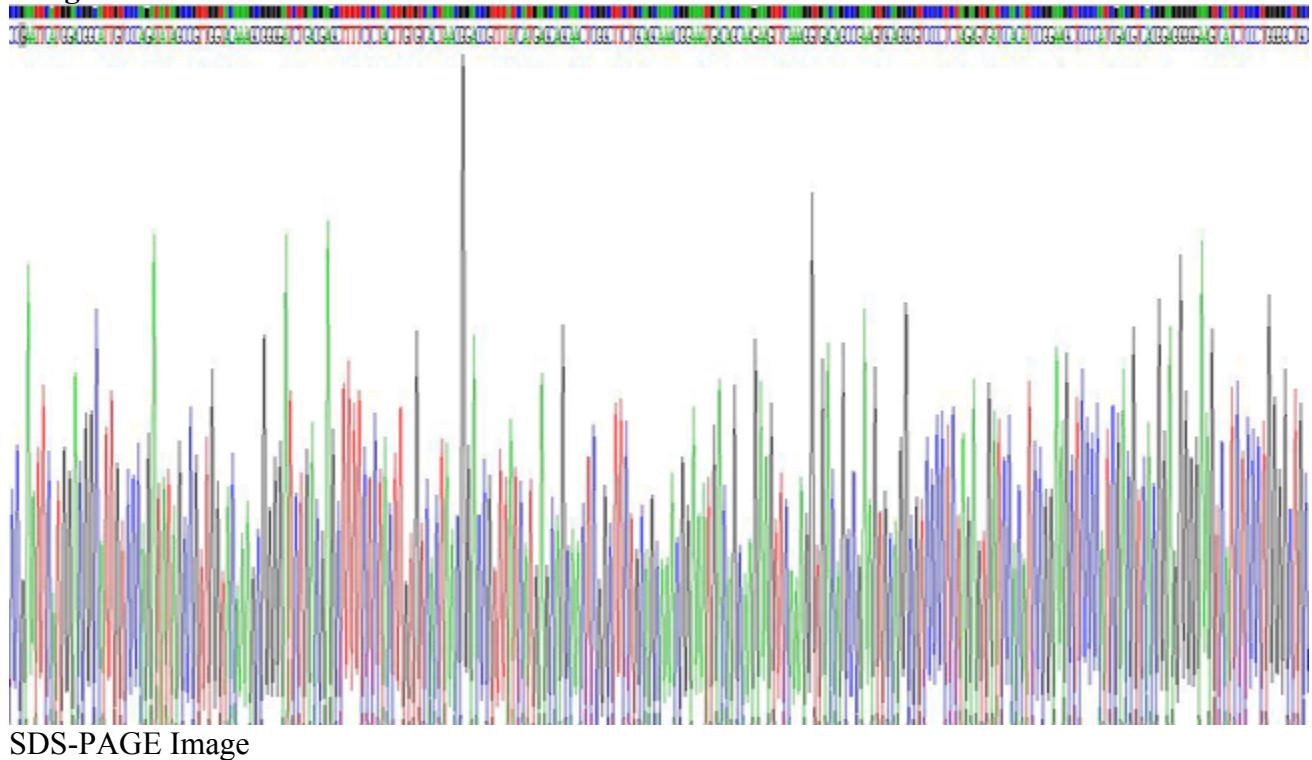

SDS-PAGE Image

### **[IMPORTANT NOTE]**

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.