Recombinant Dicer 1, Ribonuclease Type III (DICER1) Instruction Manual

SIPK128Hu01

Homo sapiens (Human)

Source Prokaryotic expression

Host E.coli

Endotoxin Level <1.0EU per 1µg (determined by the LAL method)

Subcellular LocationCytoplasmPredicted Molecular Mass35.8kDa

Accurate Molecular Mass 33kDa(Analysis of differences refer to the manual)

Residues & Tags Leu6~Ser290 with N-terminal His Tag

Buffer Formulation 100mMNaHCO3, 500mMNaCl, pH8.3, containing 0.01%

SKL, 5% Trehalose.

Traits Freeze-dried powder

Purity > 90% Isoelectric Point 5.2

Applications Positive Control; Immunogen; SDS-PAGE; WB.

SEQUENCE

LQPLS MAGLQLMTPA SSPMGPFFGL PWQQEAIHDN IYTPRKYQVE LLEAALDHNT IVCLNTGSGK TFIAVLLTKE LSYQIRGDFS RNGKRTVFLV NSANQVAQQV SAVRTHSDLK VGEYSNLEVN ASWTKERWNQ EFTKHQVLIM TCYVALNVLK NGYLSLSDIN LLVFDECHLA ILDHPYREIM KLCENCPSCP RILGLTASIL NGKCDPEELE EKIQKLEKIL KSNAETATDL VVLDRYTSQP CEIVVDCGPF TDRSGLYERL LMELEEALNF INDCNISVHS

USAGE

Reconstitute in ddH₂O to a concentration of 0.1-1.0 mg/mL. Do not vortex.


STORAGE

Avoid repeated freeze/thaw cycles. Store at 2-8°C for one month. Aliquot and store at -80°C for 12 months.

STABILITY

The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.



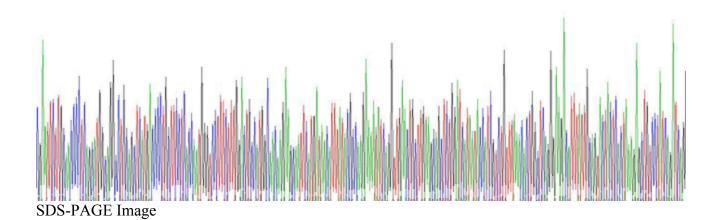

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.