

Sunlong Medical[™]

Rat TNF-α One-Step ELISA Kit

Catalog Number: One-EK015M Size: 48 Test, 96 Test

Storage: 2-8°C Validity Period: Two Years

Sensitivity:0.43 pg/mL. Assay range:4.69 pg/mL - 300 pg/mL

For the quantitative determination of rat Tumor Necrosis Factor alpha (TNF- α) concentrations in cell culture supernates, serum and plasma.

This package insert must be read entirely before using this product. For proper performance, follow the protocol provided with each individual kit.

FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.

TABLE OF CONTENTS

Introduction

Description	2
Principle of the Assay	2
Limitations of the Procedure	2

General Information

Materials Provided	3
Storage	3
Other Supplies Required	4
Precaution	4
Technical Hints	5

Assay Protocol

Sample Collection and Storage	5
Reagent Preparation	6, 7
Assay Procedure	7, 8

Analysis

Calculation of Results	
Typical Data	
Sensitivity	
Precision	
Recovery	
Linearity	
Calibration	
Sample Values	
Specificity	10

ASSAY PROCEDURE SUMMARY

- 1. Prepare all reagents and standards as directed.
- 2. Add 100 μL of 2-fold diluted *Standard* to Standard well. Add 100 μL of *Assay Bufer* (1×) to Blank well.

For Serum/Plasma samples: Add 90 μL of Assay Bufer (1×) and 10 μL sample to the sample well. For Cell culture supernates: Add 50 μL of Assay Bufer (1×) and 50 μL of

cell culture supernates to the sample well.


4. Incubate for 1 hour at 22-28°C.

5. Aspirate and wash 4 times with 60 seconds interval.

Add 100 μL of Substrate Solution to each well.
 Incubate for 10±5 minutes at 22-28 °C. Protect from light.

7. Add 100 µL of Stop Solution to each well.

8. Read at 450 nm within 30 minutes. Correction 570 or 630 nm.

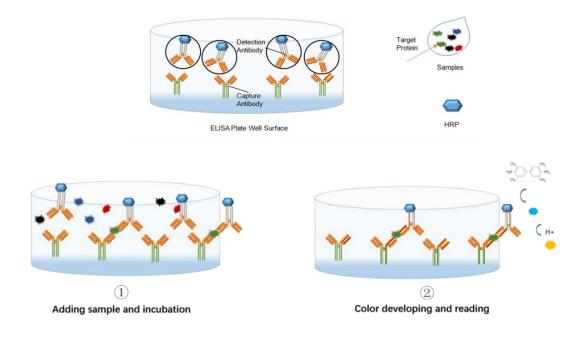



DESCRIPTION

Tumor necrosis factor alpha (TNF- α), also known as cachectin and TNFSF1A, is an adipokine involved in systemic inflammation and is a member of a group of cytokines that stimulate the acute phase reaction. It is produced chiefly by activated macrophages, although it can be produced by many other cell types such as CD4⁺ lymphocytes, NK cells, neutrophils, mast cells, eosinophils, and neurons. It functions as a multipotent modulator of immune response and further acts as a potent pyrogen. TNF- α circulates throughout the body responding to stimuli (infectious agents or tissue injury), activating neutrophils, altering the properties of vascular endothelial cells, regulating metabolic activities of other tissues, as well as exhibiting tumoricidal activity by inducing localized blood clotting. TNF- α may play a significant role in the pathogenesis of inflammatory disease of the joints and other tissues. Recently, a growing body of information has pointed to a role for TNF- α in the pathogenesis of AIDS. Measurement of TNF- α levels has also been shown to be useful in transplant research.

PRINCIPLE OF THE ASSAY

This assay employs the quantitative sandwich enzyme immunoassay technique. A monoclonal antibody specific for Rat TNF- α has been immobilized onto microwells, and two pellets of the biotin-linked detect antibody specific for TNF- α (light yellow) and streptavidin-HRP (purple) are pre-placed in the microwells, sealed by the adhesive film. Standard or samples are pipetted into the wells, then TNF- α present is bound by the immobilized antibody and detect antibody, of which the latter is conjugated with streptavidin-HRP in the incubation. After washing, substrate solution reacts with HRP and color develops in proportion to the amount of TNF- α bound by the immobilized antibody. The color development is stopped and the intensity of the color is measured by microplate reader.

LIMITATIONS OF THE PROCEDURE

- ③ FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.
- ③ Do not use expired kit or reagents.
- ③ Do not use reagents from other lots or manufacturers. Do not prepare component by yourself.
- If concentration of assayed factor in samples is higher than the highest standard, dilute the serum/plasma samples with *Assay Bufer*, dilute the cell culture supernate samples with *cell culture medium*. Reanalyze these and multiply results by the appropriate dilution factor.
- Any variation in testing personnel, sample preparation, standard dilution, pipetting technique, washing techniques, incubation time, temperature, kit age and equipment can cause variation in results.
- Bacterial or fungal contamination in either samples or reagents, or cross-contamination between reagents may cause erroneous results.
- This assay is designed to eliminate interference by factors present in biological samples. Until all factors have been tested in the ELISA immunoassay, the possibility of interference cannot be excluded.



MATERIALS PROVIDED (96 Test)

Unopened kit should be stored at 2 - 8°C.



- TNF- α Microplate (1 plate): 96-well polystyrene microplate coated with a monoclonal antibody against Rat TNF- α , two pellets of the biotin-linked detect antibody specific for TNF- α (light yellow) and streptavidin-HRP (purple) in the microwells, one piece of plate adhesive film with the crossing opening.
- TNF- α Standard (2 vials): Recombinant Rat TNF- α in a buffered protein base with preservatives; lyophilized.
- Assay Buffer (10×) (1 bottle, 5 mL): PBS with 15 % Tween-20.
- Substrate (1 bottle, 11 mL): TMB (tetramethyl-benzidine).
- Stop Solution (1 bottle, 11 mL): 0.18 M sulfuric acid.
- S Washing Buffer (20×) (1 bottle, 11 mL): PBS with 1 % Tween-20.

STORAGE

Store kit reagents between 2 and 8°C. Immediately after use remaining reagents should be returned to cold storage (2 to 8°C). Expiry of the kit and reagents is stated on labels.

Expiration date of the kit components can only be guaranteed if the components are stored properly, and if, in case of repeated use of one component, this reagent is not contaminated by the first handling.

Unopened kit		Store at 2 - 8 $^{\circ}$ C (See expiration date on the label).		
Opened/	1× Washing Buffer 1× Assay Buffer Stop Solution Substrate TMB	Up to 1 month at 2 - 8°C.		
Reconstituted Reagents	Standard	Up to 1 month at $\leq -20^{\circ}$ C in a manual defrost freezer. Discard after use.		
Microplate Wells		Up to 1 month at 2 - 8°C. Return unused strips to the foil pouch containing the desiccant pack, reseal along entire edge to maintain plate integrity.		

Provided this is within the expiration date of the kit.



OTHER SUPPLIES REQUIRED

- Solution Microplate reader capable of measuring absorbance at 450 nm, with correction wavelength set at 570 nm or 630 nm.
- Bipettes and pipette tips.
- 3 50 µL to 300 µL adjustable multichannel micropipette with disposable tips.
- Aultichannel micropipette reservoir.
- Beakers, flasks, cylinders necessary for preparation of reagents.
- **Beionized or distilled water**.
- Output Polypropylene test tubes for dilution.

PRECAUTION

- All chemicals should be considered as potentially hazardous.
- We therefore recommend that this product is handled only by those persons who have been trained in laboratory techniques and that it is used in accordance with the principles of good laboratory practice. Wear suitable protective clothing such as laboratory overalls, safety glasses and gloves.
- Care should be taken to avoid contact with skin or eyes. In the case of contact with skin or eyes wash immediately with water. See material safety data sheet(s) and/or safety statement(s) for specific advice.
- The Stop Solution provided with this kit is an acid solution. Wear eyes, hand, face, and clothing protection when using this material.
- ③ Reagents are intended for research use only and are not for use in diagnostic or therapeutic procedures.
- ③ Do not mix or substitute reagents with those from other lots or other sources.
- ③ Do not use kit reagents beyond expiration date on label.
- ⁽³⁾ Do not expose kit reagents to strong light during storage and incubation.
- ⁽³⁾ Do not eat or smoke in areas where kit reagents or samples are handled.
- Avoid contact of skin or mucous membranes with kit reagents or specimens.
- ③ Rubber or disposable latex gloves should be worn while handling kit reagents or specimens.
- Avoid contact of substrate solution with oxidizing agents and metal.
- Avoid splashing or generation of aerosols.
- In order to avoid microbial contamination or cross- contamination of reagents or specimens which may invalidate the test use disposable pipette tips and/or pipettes.
- Itse clean, dedicated reagent trays for dispensing the conjugate and substrate reagent.
- ⁽³⁾ Exposure to acid inactivates the HRP and antibody conjugate.
- ③ Glass-distilled water or deionized water must be used for reagent preparation.
- Substrate solution must be warmed to room temperature prior to use.
- Becontaminate and dispose specimens and all potentially contaminated materials as they could contain infectious agents. The preferred method of decontamination is autoclaving for a minimum of 1 hour at 121.5 ℃.
- Liquid wastes not containing acid and neutralized waste may be mixed with sodium hypochlorite in volumes such that the final mixture contains 1.0 % sodium hypochlorite. Allow 30 minutes for effective decontamination. Liquid waste containing acid must be neutralized prior to the addition of sodium hypochlorite.
- In some cases, an insoluble precipitate of stabilizing protein has been seen in the Standard Diluent. This precipitate does not interfere in any way with the performance of the test and can thus be ignored. Or remove precipitate by centrifuging at $6,000 \times g$ for 5 minutes.



TECHNICAL HINTS

- ③ When mixing or reconstituting protein solutions, always avoid foaming.
- To avoid cross-contamination, change pipette tips between additions of each standard level, between sample additions, and between reagent additions. Also, use separate reservoirs for each reagent.
- When using an automated plate washer, adding a 30 seconds soak period before washing step and/or rotating the plate between wash steps may improve assay precision.
- ③ To ensure accurate results, proper adhesion of plate sealers during incubation steps is necessary.
- Substrate Solution should remain colorless until added to the plate. Keep Substrate Solution protected from light. Substrate Solution should change from colorless to gradations of blue.
- Stop Solution should be added to the plate in the same order as the Substrate Solution.
- Solution. Wells that are green in color indicate that the Stop Solution has not mixed thoroughly with the Substrate Solution.
- It is recommended that all samples and standards be assayed in duplicate.
- Take care not to scratch the inner surface of the microwells.
- B High background maybe occur when using some cell culture medium as Assay Buffer (1×).

SAMPLE COLLECTION AND STORAGE

Cell Culture Supernates – Remove particulates by centrifugation at $300 \times g$ for 10 minutes and assay immediately or aliquot and store samples at ≤ -20 °C for later use. Avoid repeated freeze-thaw cycles.

Serum – Use a serum separator tube (SST) and allow samples to clot for 30 minutes before centrifugation for 10 minutes at 1,000 × g. Remove serum and assay immediately or aliquot and store samples at ≤ -20 °C for later use. Avoid repeated freeze-thaw cycles.

Plasma – Collect plasma using EDTA, citrate or heparin as anticoagulant. Centrifuge at 1,000 × g within 30 minutes of collection. Assay immediately or aliquot and store samples at \leq -20°C for later use. Avoid repeated freeze-thaw cycles.

Other biological samples might be suitable for use in the assay. Cell culture supernates, serum and plasma were tested with this assay. Dilution with Assay Buffer maybe needed.

Note: Samples containing a visible precipitate must be clarified prior to use in the assay. Do not use grossly hemolyzed or lipemic specimens.

Samples should be aliquoted and must be stored frozen at -20° C to avoid loss of bioactivity. If samples are to be run within 24 hours, they may be stored at 2 to 8 °C.

Avoid repeated freeze-thaw cycles. Prior to assay, the frozen sample should be brought to room temperature slowly and mixed gently.



REAGENT PREPARATION

Bring all reagents and samples to room temperature before use.

If crystals form in the Buffer Concentrates, warm and gently stir them until completely dissolved.

Washing Buffer (1.5×)

Pour 10 mL of the **Washing Buffer (20**×) into a clean 250 mL graduated cylinder. Bring to final volume of $130 \sim 135$ mL with pure or deionized water.

Mix gently to avoid foaming.

Transfer to a clean wash bottle and store at 2 to $25 \degree$ C. Washing Buffer ($1.5 \times$) is stable for 30 days.

Assay Buffer (1×)

Pour the entire contents (5 mL) of the Assay Buffer $(10 \times)$ into a clean 100 mL graduated cylinder. Bring to final volume of 50 mL with distilled water (For serum/plasma samples) or cell culture medium (For cell culture supernates). Mix gently to avoid foaming.

Store at 2 to 8° C. Assay Buffer (1×) is stable for 30 days.

Sample Dilution


If your samples have high TNF- α content, dilute serum/plasma samples with Assay Buffer (1×).

Rat TNF-a Standard

Briefly spin the vial at 6,000 rpm for 30 seconds before opening. Reconstitute **Rat TNF-** α **Standard** by addition of **distilled water** (For serum/plasma samples) or **cell culture medium** (For cell culture supernates). Reconstitution volume is stated on the label of the standard vial. Swirl or mix gently to insure complete and homogeneous solubilization (concentration of reconstituted standard = 600 pg/mL).

Allow the standard to reconstitute for 10 - 30 minutes. Mix well prior to making dilutions. Use polypropylene tubes.

Mixing concentrated Rat TNF- α standard (230 µL) with 230 µL of Assay Bufer (1 ×) creates the high standard (300 pg/mL). Pipette 230 µL of Assay Bufer (1 ×) into each tube. Use the high standard to produce a 1:1 dilution series (scheme below). Mix each tube thoroughly before the next transfer. Assay Bufer (1 ×) serves as the zero standard (0 pg/ml).

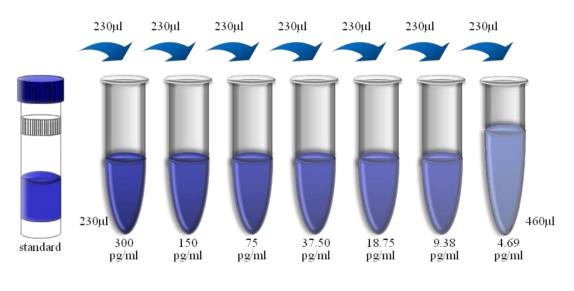



ASSAY PROCEDURE

Bring all reagents and samples to room temperature before use.

- 1. Prepare all reagents including microplate, samples, standards and working solution as described in the previous sections.
- 2. Remove excess microplate strips and return them to the foil pouch containing the desiccant pack, and reseal for further use. In any case, avoid touching the inner surface of the microwells and gently tap the plate to ensure that the pellets on the bottom of the microwells.
- 3. Adding Standard: Add 100 μ L of 2-fold diluted *Standard* to Standard well. Add 100 μ L of *Assay Bufer* (1 ×) to Blank well. *The standards/samples can be added directly from the middle opening of the plate adhesive film.*
- 4. Adding Samples: For Serum/Plasma samples Add 90 μ L of *Assay Bufer (1 ×)* and 10 μ L sample to the sample well. For Cell culture supernates Add 50 μ L of *Assay Bufer (1 ×)* and 50 μ L sample to the sample well.
- 5. Incubation: Cover with an adhesive strip. Incubate at room temperature (22-28°C) for 1 hour on a microplate shaker set at 300-500 rpm.
- 6. Washing: Removing the plate adhesive film. Aspirate each well and wash by filling each well with 300μL of *Washing Bufer (1.5 ×)*, repeating the process 3 times for a total four washes with 60 seconds interval. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Buffer by aspirating or decanting. Invert the plate and blot it against clean paper towels.
- 7. Adding Substrate Solution: Add 100 μL of *Substrate Solution* to each well. Incubate for 10±5 minutes at room temperature (22-28℃). Protect from light.
- 8. **Stopping:** Add 100 μ L of *Stop Solution* to each well. The color will turn yellow. If the color in the well is green or if the color change does not appear uniform, gently tap the plate to ensure thorough mixing.
- 9. Reading: Measure the optical density value within 30 minutes by microplate reader set to 450 nm. If wavelength correction is available, set to 570 nm or 630 nm. If wavelength correction is not available, subtract readings at 570 nm or 630 nm from the readings at 450 nm. This subtraction will correct for optical imperfections in the plate. Reading directly at 450 nm without correction may generate higher concentration than true value.



CALCULATION OF RESULTS

Average the duplicate optical density readings for each standards and sample, then subtract the average optical density value of the zero standard.

Standard Concentration as horizontal axis, optical density (OD) Value as the vertical axis, regressing the data and create a standard curve using computer software. The data may be linearized by plotting the log of the TNF- α concentrations versus the log of the OD and the best fit line can be determined by regression analysis. This procedure will produce an adequate but less precise fit of the data.


Note: The finally concentration of top standard is 300 pg/mL.

If Serum/Plasma samples have been diluted following the instruction, the final dilution factor is 10. If Cell culture supernates have been diluted following the instruction, the final dilution factor is 2. If sample have been diluted by other means, the concentration read from the standard curve must be multiplied by the appropriate dilution factor.

TYPICAL DATA

A standard curve must be run within each assay. This standard curve is provided for demonstration only.

pg/mL	0.	O.D.		Corrected
0.00	0.082	0.080	0.081	
4.69	0.114	0.122	0.118	0.037
9.38	0.158	0.168	0.163	0.082
18.75	0.210	0.231	0.221	0.140
37.50	0.334	0.322	0.328	0.247
75.00	0.596	0.570	0.583	0.502
150.00	1.000	1.130	1.065	0.984
300.00	2.162	2.052	2.107	2.026

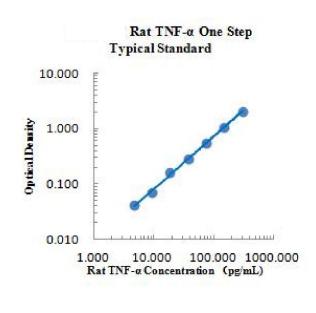



SENSITIVITY

The minimum detectable dose (MDD) of TNF- α is typically less than 0.43 pg/mL.

The MDD was determined by adding two standard deviations to the mean optical density value of ten zero standard replicates and calculating the corresponding concentration.

PRECISION

Intra-assay Precision (Precision within an assay)

Three serum-based and buffer-based samples of known concentration were tested twenty times on one plate to assess intra-assay precision.

Inter-assay Precision (Precision between assays)

Three serum-based and buffer-based samples of known concentration were tested in six separate assays to assess inter-assay precision.

	Intra-assay precision			Int	ter-assay pred	cision
Sample	1	2	3	1	2	3
n	20	20	20	6	6	6
Mean (pg/mL)	200.3	137.7	54.1	10.82	37.8	102.2
Standard deviation	17.6	13.1	4.6	0.91	2.9	6.2
CV (%)	8.8	9.5	8.4	8.4	7.7	6.1

RECOVERY

The spike recovery was evaluated by spiking 3 levels of Rat TNF- α into five health rat serum samples. The un-spiked serum was used as blank in these experiments.

The recovery ranged from 77 % to 122 % with an overall mean recovery of 95 %.

LINEARITY

To assess the linearity of the assay, five samples were spiked with high concentration of TNF- α in rat serum and diluted with Standard Diluent to produce samples with values within the dynamic range of the assay.

	Average (%)	Range (%)
1:2	99	88 - 110
1:4	104	83 - 118
1:8	105	85 - 120
1:16	93	81 - 119



CALIBRATION

This immunoassay is calibrated against a highly purified recombinant Rat TNF- α produced at SUNLONG MEDICALTM.

SAMPLE VALUES

Serum/Plasma - Thirty samples from apparently healthy rats were evaluated for the presence of $TNF-\alpha$ in this assay. No medical histories were available for the donors used in this study.

Sample Matrix	Number of Samples Evaluated	Range (pg/mL)	Detectable (%)	Mean of Detectable (pg/mL)
Serum	30	11.2-760.7	100	322.7

n.d. = non-detectable. Samples measured below the sensitivity are considered to be non-detectable.

Note: The sample range is non-physiological range. The sample range of healthy rats will difference according to geographical, ethic, sample preparation, and testing personnel, equipment varies. The above information is only reference.

SPECIFICITY

This kit could assay both natural and recombinant Rat TNF- α . A panel of substances listed below were prepared at 1 ng/mL in Standard Diluent to determine cross-reactivity. Preparations of the following substances at 1 ng/mL in a mid-range TNF- α control to determine interference. No significant cross-reactivity or interference was observed.

Human		Mouse	Rat
IL-1β	IL-18	IFN-γ	IFN-γ
IL-2	IL-21	IL-1β	IL-1β
IL-4	IL-22	IL-4	IL-4
IL-5	IL-23	IL-6	IL-6
IL-6	IFN-γ	IL-10	IL-10
IL-8	MCP-1	TNF-α	
IL-10	TGF-β1		
IL-12	TNF-α		
IL-17A	VEGF		



PLATE LAYOUT

